

 **NOTE** Inspect condition of all items before accepting shipment. File a claim with carrier if any damage is noticed.

 **CAUTION:** Before using MAXAIR Systems respirators, refer to the User's Instructions Manual received with each MAXAIR Systems Helmet.

REF 2780-06

Includes 6-Gang Battery Charger



Intended Use

The MAXAIR CAPR Cart is intended to provide storage, security, and mobile distribution of six MAXAIR Systems.

For Use With

The MAXAIR Cart is fully compatible with storage of up to six MAXAIR DLC-CAPR Systems and primary accessories.

The Cart may also be used with other MAXAIR System Configurations and select accessories.

Standard 2780-06 Cart*



*Cart is shown fully loaded. DLC-CAPR Cuff Systems and accessories are not included with the 2780-06.

Symbol Definitions



Warning, Caution, or Note



Reference Number



WARNING

Use only if package is received unopened and contents are undamaged. If damage is noted, contact the shipper for replacement or repair.

Prior to using any MAXAIR® System or component, be sure to be familiar with the system's NIOSH approved configuration.

Read and understand the User's Instructions Manual.

Failure to follow the User's Instructions Manual may be hazardous to the user's health.

The institution using this product in any application is responsible for determining the appropriateness of this equipment relative to regulatory requirements. Bio-Medical Devices Intl, Inc. does not recommend the appropriate systems for a particular institution or facility.

DO NOT use a MAXAIR System if any component is or may be damaged, contaminated, or degraded and therefore possibly unfit for safe and effective use. Replace questionable components immediately before any use of components and systems.

Use only MAXAIR Systems/ NIOSH approved compatible components.

If you need more information, contact your BMDI Sales Representative, or call BMDI customer service at 1-800-443-3842.

Materials

Cart is fabricated from powder coated galvanized steel.

Specifications

Position	Height (in/cm)	Width (in/cm)	Depth (in/cm)
Closed	43.5/110.5	36.5/92.7	21.5/54.6
Fully open	64/162.6	44/111.8	21.5/54.6
Open - front doors straight to front	64/162.6	36.5/92.7	40/101.6

Weight approximately 160 lbs./72.6 kgms (empty).



NOTE Inspect condition of all items before accepting shipment. File a claim with carrier if any damage is noticed.



WARNING

Cart weight is approximately 160 lbs. Use proper precautions and follow OSHA safety standards.



1. Safely cut and remove shipping straps.



2. Lift top outer box up and off Cart.



3. Safely cut shrink-wrapping and remove it and foam protectors from all around Cart.

4. Remove "Read-Me-First" envelope. Inside envelope are the Cart Keys and the Instructions For Use.



CAUTION

Make extra keys in the case of loss, and store in a safe place.




CAUTION

Read Instructions for Use before attempting to use the MAXAIR Cart.



5. Safely break down the sides of the bottom box.

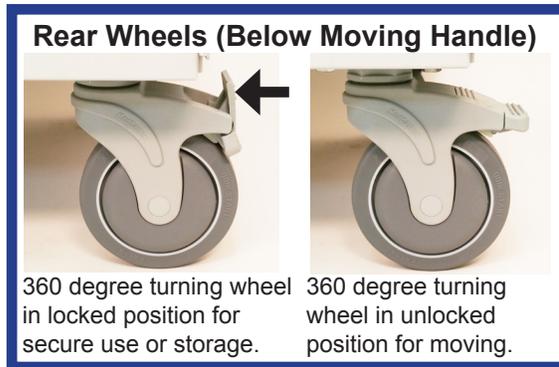

6. Unlock the rear wheels to allow the Cart to roll off the pallet. Push the locking lever forward (arrow).



7. With a slight lifting motion on the Cart Moving Handle, pull the Cart safely off the pallet and the cardboard support between the Cart and pallet.



8. Continue pulling until Cart is free from the pallet. Safely discard all wrapping materials and pallet per Institutional protocol for disposal of shipping and packing materials.